

Protocol  
K4 - Oxalic Acid Oxidation

Group F  
Jonas Adamer (12225913)  
Florian Fitsch (12218283)  
Leonhard Ritt (12208881)

Date of experiment: 2024/10/21  
Date of submission:

## Contents

<b>1</b>	<b>Objective</b>	<b>2</b>
<b>2</b>	<b>Experiment</b>	<b>2</b>
2.1	Determination of Absorbance Maximum of $\text{KMnO}_4$ . . . . .	2
2.2	Creation of Calibration Curve for $\text{KMnO}_4$ Concentration . . . . .	2
2.3	Measurement of the Autocatalyzed Redox Reaction of $\text{KMnO}_4$ and Oxalic Acid . . . .	2
2.4	Measurement of the $\text{Mn}^{2+}$ -Catalyzed Redox Reaction of $\text{KMnO}_4$ and Oxalic Acid . .	3
2.5	Measurement of the Deprotonation of Phenolphthalein by Sodium Hydroxide . . . . .	3
<b>3</b>	<b>Results and Discussion</b>	<b>4</b>
3.1	Absorbance Maximum and Calibration Curve of $\text{KMnO}_4$ . . . . .	4
3.2	Reaction order of Autocatalyzed Redox Reaction of $\text{KMnO}_4$ and Oxalic Acid . . . . .	5

## 1 Objective

In this assignment, the kinetics of chemical reactions are studied through two separate experiments.

In the first experiment, the redox reaction between oxalic acid and potassium permanganate is studied: After determining the absorbance maximum of potassium permanganate, a series of solutions is produced and measured as a calibration for concentration of permanganate as a function of absorbance. Finally, the change in permanganate-concentration is measured in-situ, both with and without a catalyst, in order to determine the reaction order.

In the second experiment, the acid-base reaction between phenolphthalein and sodium hydroxide is studied. The absorbance of the reaction mixture is also measured in-situ and the reaction order regarding phenolphthalein determined.

## 2 Experiment

### 2.1 Determination of Absorbance Maximum of $\text{KMnO}_4$

In order to find the absorbance maximum, 79 mg of potassium permanganate were weighed out, dissolved in deionized water and diluted in a 50 mL volumetric flask, yielding a 0.01 mol/L solution. 5 mL of this solution were taken out and once again diluted in a 50 mL flask, yielding a concentration of 1 mmol/L.

2 mL of the prepared solution were filled into a cuvette and inserted into the photometer. Its absorbance was then measured in wavelength increments of 10 nm between 450 nm and 600 nm. For each wavelength, the photometer was calibrated by setting values of 0% and 100% transmission, using an opaque block and a cuvette filled with deionized water, respectively. Since high absorbances were measured at 530 nm and 550 nm, another series of measurements was taken between 520 nm and 560 nm with increments of 2 nm. Using this process, a maximum at 532 nm was determined.

### 2.2 Creation of Calibration Curve for $\text{KMnO}_4$ Concentration

To create the calibration, the 1 mmol/L potassium permanganate solution, which was prepared in Section 2.1, was further diluted to prepare six further solutions, with concentrations of 0.75, 0.60, 0.50, 0.40, 0.25 and 0.10 mmol/L. After calibrating the photometer as described in Section 2.1, the absorbance values of the seven solutions were measured at the previously determined wavelength of 532 nm. To account for inaccuracies, each measurement was taken three times, with the cuvette being filled with fresh solution after every measurement.

### 2.3 Measurement of the Autocatalyzed Redox Reaction of $\text{KMnO}_4$ and Oxalic Acid

By further diluting the 10 mmol/L solution prepared in 2.1, a solution with a concentration of 4 mmol/L was prepared. Additionally, a 40 mmol/L solution of oxalic acid was obtained by dissolving 4.000 g oxalic acid in deionized water in a 100 mL volumetric flask. Additionally, a sulfuric acid solution of about 24 w% was prepared by combining three parts water with one part concentrated (96 %) sulfuric acid.

To start the reaction, 10 mL of the oxalic acid solution, 0.4 mL of the sulfuric acid and 4 mL of the potassium permanganate solution were added to a beaker and mixed thoroughly. A stopwatch was started once the permanganate solution was added. After filling 2 mL of the mixture into a cuvette, absorbance measurements were taken every 10 seconds, the measured value reached 0 or stopped changing. This experiment was repeated two more times.

## 2.4 Measurement of the $\text{Mn}^{2+}$ -Catalyzed Redox Reaction of $\text{KMnO}_4$ and Oxalic Acid

By dissolving AAAAAA g manganese sulfate in deionized water in a AAAAAA mL volumetric flask, a 40 mmol/L solution was prepared.

Three reactions were measured as described in Section 2.3, with the addition of adding 0.2 mL of the prepared manganese sulfate solution to the reaction mixture before measuring it.

## 2.5 Measurement of the Deprotonation of Phenolphthalein by Sodium Hydroxide

By diluting 7.5 mL of 2 mol/L sodium hydroxide with deionized water in a 50 mL volumetric flask, 50 mL of a 0.3 mol/L solution was prepared. Additionally, by dissolving 877 g of sodium chloride in deionized water and diluting to 50 mL, a 0.3 mol/L solution was prepared.

By diluting the former with the latter solution, small amounts of solutions with sodium hydroxide concentrations of 0.2, 0.1 and 0.05 mol/L, were prepared.

2 mL of each of the four solutions were then transferred into a cuvette, and a drop of phenolphthalein was added. The cuvettes were inserted into the photometer, which was set at 532 nm<sup>1</sup>, and, after dropping to an absorbance value of 0.8, further measurements were taken every 30 seconds over a span of six minutes. Since the time available for the experiment ran out, only one sample of each concentration was measured. The final sample (0.05 mol/L) did not reach the starting absorbance of 0.8 in time, and as such, higher absorbance values were recorded and used for the evaluation of this sample.

---

<sup>1</sup>Due to an error, the desired wavelength of 550 nm was not set before starting the experiment

### 3 Results and Discussion

#### 3.1 Absorbance Maximum and Calibration Curve of $\text{KMnO}_4$

The recorded absorbances of potassium permanganate are plotted as an absorbance spectrum in Figure 1. The determined absorbance maximum is at 532 nm.

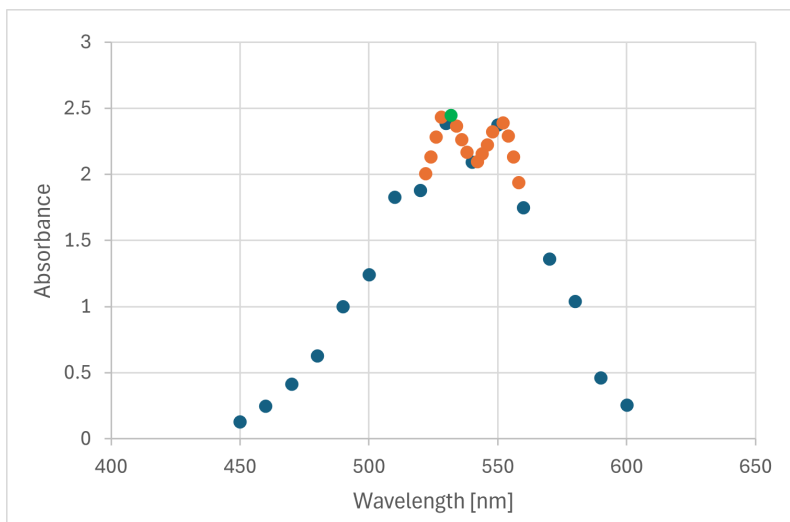


Figure 1: Absorbance spectrum of  $\text{KMnO}_4$  at wavelengths between 450 and 600 nm. The blue dots show the measurements of the first scan, which was taken in 10 nm increments, while the orange dots mark the measurements of the second scan (2 nm increments). The absorbance maximum at 532 nm is shown in green.

Figure 2 shows the calibration curve (created through linear regression) for the concentration of potassium permanganate as a function of absorbance. As can be seen, linearity, and thus the validity of the Lambert-Beer law is given for the entire tested concentration range.

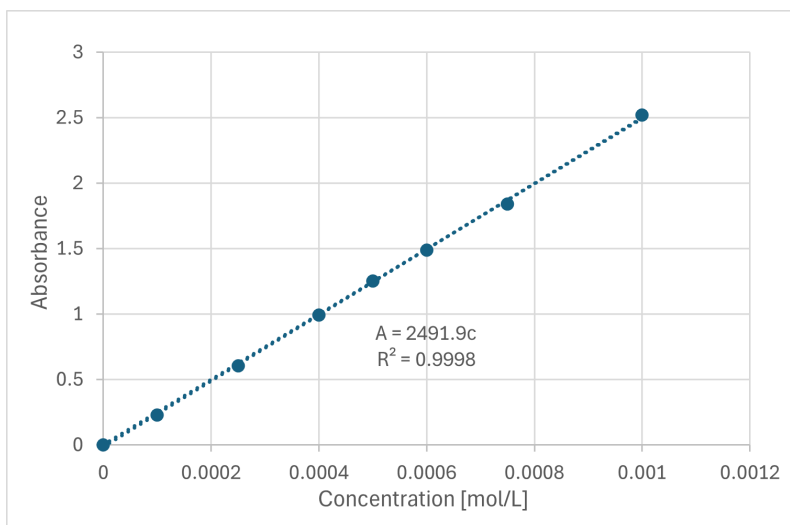


Figure 2: Average values of measured absorbance plotted against potassium permanganate concentration. A calibration curve is fit to the data using linear regression.

By inserting the resulting linear equation  $A = 2491.9 \cdot c$  into the Beer-Lambert law (Equation 1), rearranging and inserting the known cuvette length of 1 cm, the molar attenuation coefficient  $\varepsilon$  can be calculated:

$$A = \varepsilon \cdot c \cdot d = 2491.9 \cdot c \quad (1)$$

$$\varepsilon = \frac{2491.9}{d} = 2491.9 \quad (2)$$

where:  $A$  = Absorbance

$\varepsilon$  = Molar attenuation coefficient [ $\text{L mol}^{-1} \text{ cm}^{-1}$ ]

$c$  = Concentration of sample [ $\text{mol/L}$ ]

$d$  = Path length through sample (cuvette length) [cm]

### 3.2 Reaction order of Autocatalyzed Redox Reaction of $\text{KMnO}_4$ and Oxalic Acid

By rearranging the Lambert-Beer law (Equation 1), the concentration of a sample can be determined from a given absorption value.

$$c = \frac{A}{\varepsilon \cdot d} = \frac{A}{2491.9} \quad (3)$$